

Electrocardiograph



E3A

3-Channel Digital Electrocardiograph

Features


- 7" high resolution TFT LCD display
- 12 leads simultaneous acquisition
- One-touch operation with user definable interface
- Full alphanumeric keyboard with direct function keys
- Brand new distinct reading experience with 110mm roll & Z-fold thermal paper, optional for 80mm
- Up to 3000 patient records for internal data storage, support USB, SD card, Mobile HDD
- High frequency response up to 250HZ which is exclusive module for children and neonate analysis
- Built-in rechargeable Lithium battery which can support for 2.5h continuous working

> New Elegant Design:

- 7 inches high resolution LCD display, easy to observe at each angle
- Full alphanumeric keyboard for waterproof
- Direct fast keys for easy operation
- Built-in rechargeable lithium battery for more than 2.5 hours continuous working
- 110mm size Roll and Z-fold thermal paper for printing
- Multi ports: USB, SD card, LAN



Waterproof keyboards

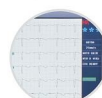

Color display



Safety package



Multi ports

> Large Data Manage:

- Internal memory reaches to 3000 patient files.
- Optional for USB, SD card, Mobile HDD for more memory storage
- ECG management software for optional

> Reliable Analysis

- Advanced interpretation system with CEIS database
- New platform application
- 12 leads simultaneous acquisition

> High frequency response

- 0.05Hz-250Hz frequency response
- 0.05Hz for high-pass frequency, which is more sensitive for weaker ECG signal detection.
- 250Hz for low-pass frequency, which is special for neonate and children Smart workflow
- Sound and Visual alarm
- Multi waveform display mode
- Working mode including automatic, manual, rhythm
- Multi format: XML, JPG, Dicom
- Multi language for different choice

Somo International Co., Ltd

Shenzhen Office: Floor 7, Jinyuan Commercial Building, No. 300 Xixiang Dadao, Baoan District, Shenzhen, China
Tel : 0086-755 2772 8713 Fax : 0086-755 2772 5003 E-mail : sales@somo-med.com info@somo-med.com Website : www.somo-med.com

Electrocardiograph



E6

6- Channel Digital Electrocardiograph

Features

- Compact and portable design, light in weight
- 7-inch high resolution LCD with touch screen, easy operation
- Waterproof full alphanumeric silicon keyboard
- 12 leads simultaneous acquisition
- Brand new distinct reading experience with 110mm roll & Z-fold thermal paper
- Massive data storage capacity which can up to 3000 patient records for internal data storage
- Built-in rechargeable Lithium battery which can support for 2.5h continuous working
- Wireless connection for optional
- Various data formats: JPG, XML, DICOM and so on

> New Elegant Design:

- 7 inches high resolution LCD display, easy to observe at each angle
- One-touch operation Touch screen
- Full alphanumeric keyboard for waterproof
- Direct fast keys for easy operation
- Built-in rechargeable lithium battery for more than 2.5 hours continuous working
- 110mm size Roll and Z-fold thermal paper for printing
- Multi ports: USB, SD card

> Large Data Manage:

- Internal memory reaches to 3000 patient files.
- Optional for USB, SD card, Mobile HDD for more memory storage
- ECG management software for optional

> Reliable Analysis

- Advanced interpretation system with CEIS database
- New platform application
- 12 leads simultaneous acquisition
- Wifi connection with ECG management software

> High frequency response

- 0.05Hz-250Hz frequency response
- 0.05Hz for high-pass frequency, which is more sensitive for weaker ECG signal detection.
- 250Hz for low-pass frequency, which is special for neonate and children

> Smart workflow

- Sound and Visual alarm
- Multi waveform display mode
- Working mode including automatic, manual, rhythm
- Multi format: XML, JPG, Dicom
- Multi language for different choice



Color display screen



Waterproof keyboards



Ports



Portable Carrying Handle



Elegant design

Somo International Co., Ltd

Shenzhen Office: Floor 7, Jinyuan Commercial Building, No. 300 Xixiang Dadao, Baoan District, Shenzhen, China
Tel : 0086-755 2772 8713 Fax : 0086-755 2772 5003 E-mail : sales@somo-med.com info@somo-med.com Website : www.somo-med.com

Electrocardiograph



E12A

12- Channel Digital Electrocardiograph

Features

- ▶ Compact and portable design, light in weight
- ▶ 10.4 inch high resolution LCD with touch screen, easy operation
- ▶ Waterproof full alphanumeric silicon keyboard
- ▶ 12 leads simultaneous acquisition
- ▶ A4 size 210mm Z-fold thermal paper
- ▶ Massive data storage capacity which can up to 3000 patient records for internal data storage
- ▶ Built-in rechargeable Lithium battery which can support for 2.5h continuous working
- ▶ Wireless connection for optional
- ▶ Various data formats: JPG, XML, DICOM and so on



Color display screen



Waterproof keyboards

New Elegant Design:

- ▶ 10.4 inches high resolution LCD display, easy to observe at each angle
- ▶ One-touch operation Touch screen
- ▶ Full alphanumeric keyboard for waterproof
- ▶ Direct fast keys for easy operation
- ▶ Built-in rechargeable lithium battery for more than 2.5 hours continuous working
- ▶ 210mm A4 size Z-fold thermal paper for printing
- ▶ Multi ports: USB, SD card

Large Data Manage:

- ▶ Internal memory reaches to 3000 patient files.
- ▶ Optional for USB, SD card, Mobile HDD for more memory storage
- ▶ ECG management software for optional

Reliable Analysis

- ▶ Advanced interpretation system with CEIS database
- ▶ New platform application
- ▶ 12 leads simultaneous acquisition
- ▶ Wifi connection with ECG management software



Ports

High frequency response

- ▶ 0.05Hz-250Hz frequency response
- ▶ 0.05Hz for high-pass frequency, which is more sensitive for weaker ECG signal detection.
- ▶ 250Hz for low-pass frequency, which is special for neonate and children

Smart workflow

- ▶ Sound and Visual alarm
- ▶ Multi waveform display mode
- ▶ Working mode including automatic, manual, rhythm
- ▶ Multi format: XML, JPG, Dicom
- ▶ Multi language for different choice



Portable Carrying Handle

Somo International Co., Ltd

Shenzhen Office: Floor 7, Jinyuan Commercial Building, No. 300 Xixiang Dadao, Baoan District, Shenzhen, China

Tel : 0086-755 2772 8713 Fax : 0086-755 2772 5003 E-mail : sales@somo-med.com info@somo-med.com Website : www.somo-med.com